


Catheter Dressing

IV3000 Application

IV3000 10cm x 12cm

Suitable for dressing central line catheters

Step A

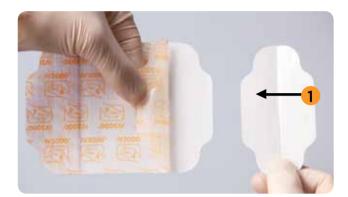

Clean and thoroughly dry the application area.

First remove section 1 from the end of the dressing. This contains 2 strips and a date label.

Step D



Using the tab gently peel off the top carrier marked 4.

Step B



The strips can be used on the catheter wings or body for improved security, to secure tubing or to stabilise catheter lumens.

Peel off the backing paper marked 2.

Step E



Mould the film around the catheter, ensuring the dressing is secure and the patient is comfortable.

Finally place the date label from section 1 on the skin or the dressing (without obscuring the view of the catheter site).

Step C



Centre the dressing over the catheter insertion site.

Place the dressing on the catheter hub and adhere the dressing to the patients skin.

Remove the remaining protector paper marked 3.



Removal

Take one corner of the film and gently pull parallel to the skin (stabilising the patient's skin and catheter with the other hand).

Release and repeat the process.

Ordering codes

Dressing						
Size	5cm x 6cm	6cm x 7cm	7cm x 9cm	9cm x 12cm	10cm x 12cm	11cm x 14cm
S&N Code	66004011	4007	4006	66004009	4008	66800512
NHS Code	ELW011	ELW032	ELW046	ELW112	ELW054	ELW363
Carton	100	100	100	50	50	25

Wound Management Smith & Nephew Healthcare Ltd Healthcare House 101 Hessle Road Hull HU3 2BN

T 01482 222200 F 01482 222211 www.smith-nephew.com/wound www.iv3000.com


IV3000 Application

IV3000 Ported 11cm x 14cm

Suitable for PICC and midline catheters


Application with a stabilisation device

Step A

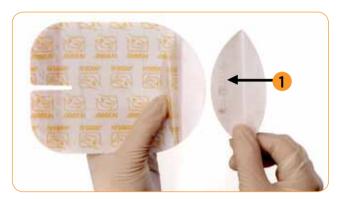

Clean and thoroughly dry the application area. First remove section 1 from the rounded end of the dressing. This contains 2 strips and a date label.

Step B

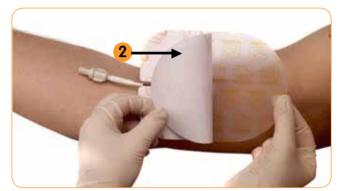


Peel off the backing paper marked 2

Step C



Place the dressing over the cannula site, ensuring that the entry site is as close to the centre of the dressing (clear section) as possible.

Step D

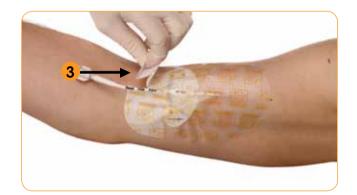

Arrange the ported area around the lumen and remove the backing papers marked 3 from the non-woven wings.

Step E



Using the tab gently peel off the top carrier marked 4.

Mould the film around the catheter, ensuring that the dresing is secure and the patient is comfortable.

Step F



The securing strips from section 1 can be used on the catheter for improved security to secure tubing or to stabilise catheter lumens.

Finally place the date label from section 1 on the skin or the dressing (without obscuring the view of the catheter site).



Removal

Gently peel back the non-woven wings to the point at which they reach the film section (stabilising the patient's skin and catheter with the other hand).

Take one corner of the film and gently pull parallel to the skin.

Release and repeat the process.

Ordering codes

Dressing						
Size	5cm x 6cm	6cm x 7cm	7cm x 9cm	9cm x 12cm	10cm x 12cm	11cm x 14cm
S&N Code	66004011	4007	4006	66004009	4008	66800512
NHS Code	ELW011	ELW032	ELW046	ELW112	ELW054	ELW363
Carton	100	100	100	50	50	25

Wound Management Smith & Nephew Healthcare Ltd Healthcare House 101 Hessle Road Hull HU3 2BN

T 01482 222200 F 01482 222211 www.smith-nephew.com/wound www.iv3000.com




IV3000 Application

IV3000 Ported 5cm x 6cm

Suitable for dressing paediatric catheters

Step A

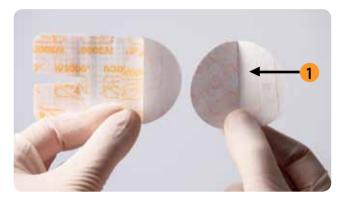

Clean and thoroughly dry the application area.

First remove section 1 from the rounded end of the dressing.

This contains 2 strips and a date label.

Step D



Position the dressing over the catheter port. Remove the backing papers marked 3 from the non woven sections and adhere over the catheter wings.

Using the tab gently peel off the top carrier marked 4.

Step B



The strips can be used on the catheter wings or body for improved security, to secure tubing or to stablise catheter lumens.

Step C



Remove the backing paper marked 2.

Step E



Mould the film around the catheter, ensuring that the dressing is secure and the patient is comfortable.

1 on the skin or the dressing (without obscuring the view of the catheter site).



Removal

Gently peel back the non-woven wings to the point at which they reach the film section (stabilising the patient's skin and catheter with the other hand).

Take one corner of the film and gently pull parallel to the skin.

Release and repeat the process.

Ordering codes

Dressing						
Size	5cm x 6cm	6cm x 7cm	7cm x 9cm	9cm x 12cm	10cm x 12cm	11cm x 14cm
S&N Code	66004011	4007	4006	66004009	4008	66800512
NHS Code	ELW011	ELW032	ELW046	ELW112	ELW054	ELW363
Carton	100	100	100	50	50	25

Wound Management Smith & Nephew Healthcare Ltd Healthcare House 101 Hessle Road Hull HU3 2BN

T 01482 222200 F 01482 222211 www.smith-nephew.com/wound www.iv3000.com


IV3000 Application

IV3000 6cm x 7cm

Suitable for dressing peripheral catheters




Step A

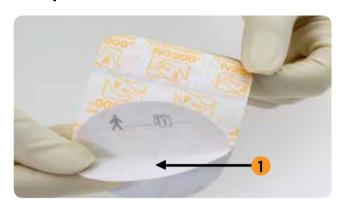

Clean and thoroughly dry the application area. First remove section 1 from the rounded end of the dressing.

This contains 2 strips and a date label.

Step D

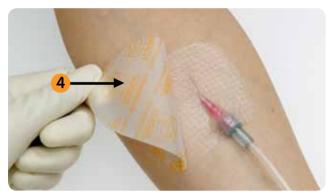


Using the tab gently peel off the top carrier marked 4.

Mould the film around the catheter to complete dressing application.

Finally place the date label from section 1 on the skin or the dressing (without obscuring the view of the catheter site).


Step B

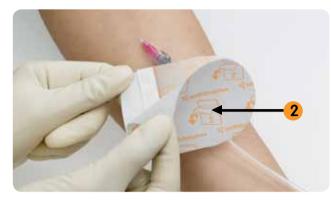

The strips can be used on the catheter wings or body for improved security, to secure tubing or to stablise catheter lumens.

Remove the backing paper marked 2.

Step C

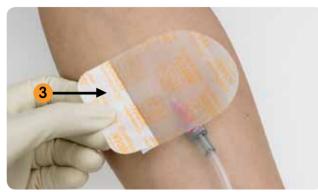

Centre the dressing over the catheter insertion

Place the dressing on the catheter port and adhere the dressing to the patients skin.

Remove the remaining protector paper marked 3.





Removal

Take one corner of the film and gently pull parallel to the skin. Release and repeat the process.

Ordering codes

Dressing						
Size	5cm x 6cm	6cm x 7cm	7cm x 9cm	9cm x 12cm	10cm x 12cm	11cm x 14cm
S&N Code	66004011	4007	4006	66004009	4008	66800512
NHS Code	ELW011	ELW032	ELW046	ELW112	ELW054	ELW363
Carton	100	100	100	50	50	25

Wound Management Smith & Nephew Healthcare Ltd Healthcare House 101 Hessle Road Hull HU3 2BN

T 01482 222200 F 01482 222211

www.smith-nephew.com/wound www.iv3000.com



IV3000 Application

IV3000 Ported 7cm x 9cm

Suitable for dressing peripheral catheters

Step A



Clean and thoroughly dry the application area. First remove section 1 from the rounded end of the dressing.

This contains 2 strips and a date label.

Step B



The strips can be used on the catheter wings or body for improved security, to secure tubing or to stablise catheter lumens.


Remove the backing paper marked 2.

Step C



Position the dressing over the catheter port. Remove the backing papers marked **3** from the non woven sections and adhere over the catheter wings.

Step D

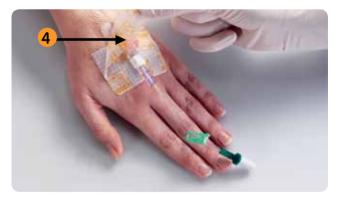

Using the tab, gently peel off the top carrier marked 4

Step E



Mould the film around the catheter, ensuring that the dressing is secure and the patient is comfortable.

Finally place the date label from section 1 on the skin or the dressing (without obscuring the view of the catheter site).

Removal

Gently peel back the non-woven wings to the point at which they reach the film section (stabilising the patient's skin and catheter with the other hand).

Take one corner of the film and gently pull parallel to the skin. Release and repeat the process.

Ordering codes

Dressing						
Size	5cm x 6cm	6cm x 7cm	7cm x 9cm	9cm x 12cm	10cm x 12cm	11cm x 14cm
S&N Code	66004011	4007	4006	66004009	4008	66800512
NHS Code	ELW011	ELW032	ELW046	ELW112	ELW054	ELW363
Carton	100	100	100	50	50	25

Wound Management Smith & Nephew Medical Ltd 101 Hessle Road Hull HU3 2BN

T +44 (0) 1482 225181 F +44 (0) 1482 328326 www.smith-nephew.com/wound www.iv3000.com




Catheter Dressing

IV3000 Application

IV3000 Ported 9cm x 12cm

Suitable for dressing jugular and central line catheters


Step A

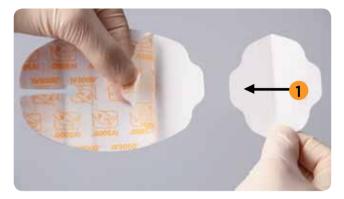

Clean and thoroughly dry the application area. First remove section 1 from the rounded end of the dressing.

This contains 2 strips and a date label.

Step D

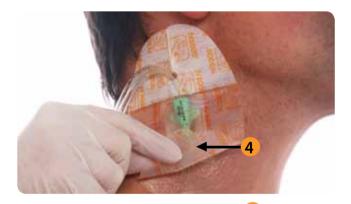

Peel off the top film carrier marked 4 using the tab at the base of the dressing, pulling gently towards the direction of the port.

Step B



The strips can be used on the catheter wings or body for improved security, to secure tubing or to stabilise catheter lumens.

Peel off the backing paper marked **2** Place the dressing over the cannula site, ensuring that the entry site is as close to the centre of the dressing (clear section) as possible.

Step E



Mould the film around the catheter, ensuring that the dressing is secure and the patient is comfortable.

Finally place the date label from section 1 on the skin or the dressing (without obscuring the view of the catheter site).

Step C



Arrange the ported area around the lumen and remove the backing papers marked 3 from the non-woven wings.

The non-woven wings can be gently pulled together and overlapped to ensure the weight of the lumen is well secured.




Removal

Gently peel back the non-woven wings to the point at which they reach the film section (stabilising the patient's skin and catheter with the other hand).

Then take one corner of the film and gently pull parallel to the skin.

Release and repeat process.

Ordering codes

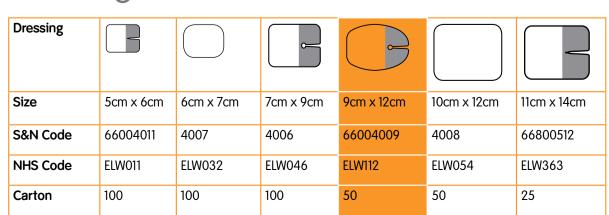

Wound Management Smith & Nephew Healthcare Ltd Healthcare House 101 Hessle Road Hull HU3 2BN

T 01482 222200 F 01482 222211 www.smith-nephew.com/wound www.iv3000.com