

RUSCH CARE

URINE COLLECTORS



Teleflex
MEDICAL

SIMPLE, SAFE AND HYGIENIC:

With their economic and innovative high-quality products, Teleflex Medical are improving the quality of life of their urinary drainage patients and are offering them maximum protection against nosocomial urinary tract infections. About 90% of these infections are caused by conventional urinary drainage. Our closed urinary drainage system, which meets all the hygienic demands made on long-term medical care, represents the first and most important prophylactic measure.



**URINE SAMPLING
WITHOUT THE NEED
FOR A NEEDLE**

Urine sampling is carried out at the ridged universal connector with a Luer syringe without a needle. Thus, needle stick injuries are eliminated and the risk of invading bacteria is greatly reduced.

S-BAG - BED URINE COLLECTOR

**2000 ML, CLOSED URINARY
DRAINAGE SYSTEM FOR
STANDARD DRAINAGE
UNDER LONG-TERM
CONDITIONS**

**Hygienic and quick emptying**

Non-dripping bottom outlet with one-hand clamp.
Doghouse holder to fix the outlet tube, to prevent contamination and to achieve additional floor clearance.

REF.	VOLUME	TUBE	QTY	BOX
85 04 14	2000 ml	120 cm	1	20

- connector with needle-free urine sampling port, with cap
- drip chamber with integrated check valve
- double bacteria-proof and liquid-tight ventilation
- non-dripping bottom outlet with one-hand clamp
- doghouse holder to fix the outlet tube
- clear graduation
- flexible and non-kinking tube, with sliding clamp and sheet clip, extra cord with loop
- latex-free
- double packed
- **sterile, disposable**

**Vented drip chamber with integrated check valve**

The vented drip chamber with high drip level and safely functioning valve reliably prevents urine reflux and intraluminally ascending infections. Optimal pressure compensation and safe outlet performance thanks to bacteria-proof ventilation. Safe even under high pressure load.



Safe centred fixation
with double hook.
Additional safety thanks
to extra cord with loop.

BED URINE COLLECTORS

2000 ML

REF.	VOLUME	TUBE	QTY	BOX
85 04 12	2000 ml	100 cm	75	75
85 04 19	2000 ml	120 cm	75	75

- connector with needle-free urine sampling port and with cap
- with non-reflux valve
- with T-tap outlet
- clear graduation
- flexible and non-kinking tube, with clamp
- without hanger
- latex-free
- **sterile, disposable**



BED URINE COLLECTORS

2000 ML

REF.	VOLUME	TUBE	QTY	BOX
85 04 07	2000 ml	90 cm	75	75
85 04 18	2000 ml	130 cm	75	75

- connector with cap
- with non-reflux valve
- with T-tap outlet
- clear graduation
- flexible and non-kinking tube
- without hanger
- latex-free
- **sterile, disposable**



URINE SAMPLING WITHOUT THE NEED FOR A NEEDLE

Urine sampling is carried out at the ridged universal connector with a Luer syringe without a needle. Thus, needle-stick injuries are eliminated and the risk of invading bacteria is greatly reduced.



T-TAP OUTLET

Easy to use; it can even be operated using one hand only.

HANGER

Firm fixation at bedside.



HANGER			
alone, for bed urine collectors			
REF.	QTY	BOX	
85 04 45	10	180	



CAUTION:

Non-sterile urine collectors should not be connected with sterile trans-urethral, suprapubic or nephrostomy catheters or any other sterile catheter. Non-sterile urine collectors can be connected only with non-sterile drainage devices (e.g. external catheters).

BED URINE COLLECTORS**2000 ML**

REF.	VOLUME	TUBE	QTY	BOX
85 04 16	2000 ml	90 cm	30	240
85 04 00	2000 ml	90 cm	250	250
85 04 17*	2000 ml	90 cm	30	210
85 04 15	2000 ml	130 cm	30	210
85 04 02	2000 ml	130 cm	200	200
85 03 46*	2000 ml	130 cm	30	210

- connector with cap
 - without outlet
 - clear graduation
 - flexible and non-kinking tube
 - with 1 hanger/Qty
 - latex-free
 - **non-sterile, disposable**
- * with non-reflux valve

**BED URINE COLLECTORS****2000 ML**

REF.	VOLUME	TUBE	QTY	BOX
85 04 01	2000 ml	90 cm	30	240
85 04 08	2000 ml	90 cm	200	200
85 04 10	2000 ml	130 cm	150	150

- connector with cap
- with non-reflux valve
- with T-tap outlet
- clear graduation
- flexible and non-kinking tube
- with 1 hanger/Qty
- latex-free
- **non-sterile, disposable**

**T-TAP OUTLET**

Easy to use; it can even be operated using one hand only.

HANGER

alone, for bed urine collectors


REF.	QTY	BOX
85 04 45	10	180





3-CHAMBER LEG URINE COLLECTORS

500 ML, WITH SKIN-FRIENDLY REVERSE SIDE

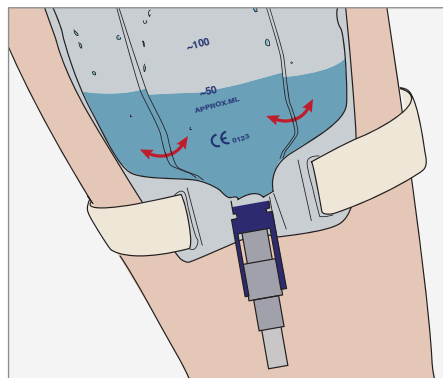

REF.	VOLUME	TUBE	QTY/BOX
85 03 50	500 ml	30 cm	8
85 03 51	500 ml	10 cm	8
85 03 52	500 ml	45 cm	8

The urine collector's division into three chambers, which are connected with each other, adds to the patient's comfort as it prevents the urine from sloshing around – also when walking. The collector fills up noiselessly and in a less visible fashion.

- connector with needle-free urine sampling port, with cap
- with non-reflux valve
- with Flip Flow outlet
- clear graduation
- flexible and non-kinking tube
- with 2 Velcro straps/Box
- latex-free
- **sterile, disposable (85 03 50, 85 03 51)**
- **non-sterile, disposable (85 03 52)**

Each pack containing a 3-chamber leg urine collector also includes two Velcro silicone-covered straps to ensure the urine collector is securely affixed.




VELCRO STRAPS

1 pair, on their own, for leg urine collectors

REF.	QTY
85 03 60 - 42 00 25	2

CAUTION:

Non-sterile urine collectors should not be connected with sterile trans-urethral, suprapubic or nephrostomy catheters or any other sterile catheter. Non-sterile urine collectors can be connected only with non-sterile drainage devices (e.g. external catheters).

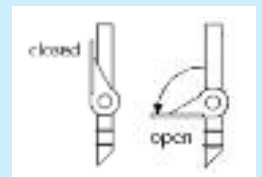

FLIP FLOP OUTLET

Extremely safe and easy to use; it can even be operated using one hand only.

A universal connector is fitted to the end of the valve, allowing for its easy connection to subsequent urine collectors.



CAUTION:

Non-sterile urine collectors should not be connected with sterile trans-urethral, suprapubic or nephrostomy catheters or any other sterile catheter. Non-sterile urine collectors can be connected only with non-sterile drainage devices (e.g. external catheters).



STANDARD STRAPS



LEG URINE COLLECTORS

750 ML

REF.	VOLUME	TUBE	QTY	BOX
85 03 30	750 ml	10 cm	10	240
85 03 31	750 ml	10 cm	30	210
85 03 32	750 ml	30 cm	10	240
85 03 33	750 ml	30 cm	30	210

- connector with cap
- with non-reflux valve
- push-pull outlet
- clear graduation
- flexible and non-kinking tube
- with 2 standard straps/Qty
- collector latex-free
- **non-sterile, disposable**

LEG URINE COLLECTORS, 750 ML

WITH SKIN-FRIENDLY REVERSE SIDE

REF.	VOLUME	TUBE	QTY	BOX
85 03 42	750 ml	10 cm	30	210
85 03 43	750 ml	30 cm	30	210

- connector with cap
- with non-reflux valve
- push-pull outlet
- clear graduation
- flexible and non-kinking tube
- with 2 standard straps/Qty
- collector latex-free
- **non-sterile, disposable**



VELCRO STRAPS

1 pair, on their own, for leg urine collectors

REF.	QTY
85 03 60 - 42 00 25	2

LEG URINE COLLECTORS

750 ML, WITH SKIN-FRIENDLY REVERSE SIDE

REF.	VOLUME	TUBE	QTY/BOX
85 03 39	750 ml	10 cm	10
85 03 41	750 ml	30 cm	10

- connector with needle-free urine sampling port, with cap
- with non-reflux valve
- with Flip Flow outlet
- clear graduation
- flexible and non-kinking tube
- with 2 Velcro straps/Qty
- collector latex-free
- **sterile, disposable**



URINE SAMPLING WITHOUT THE NEED FOR A NEEDLE

Urine sampling is carried out at the ridged universal connector with a Luer syringe without a needle. Thus, needle-stick injuries are eliminated and the risk of invading bacteria is greatly reduced.

BELLY BAG

1000 ML, FOR ALL MOBILE PATIENTS WHO ARE DISCHARGED WITH A BLADDER CATHETER

REF.	VOLUME	TUBE	QTY/BOX
B 1000	1000 ml		10
B 1000 P	1000 ml		10
	with needle-free urine sampling port		
B 1000 CT	1000 ml	60 cm	10

REPLACEMENT BELT	
REF.	QTY/BOX
PO 337	10

- skin-friendly and soft reverse side
- non-dripping bottom outlet with twist-off tap enables hygienic and quick emptying of the bag
- with non-reflux valve
- handy retaining band
- Ref. B 1000 P: connector with needle-free urine sampling port
- clear graduation
- latex-free
- **sterile, disposable**



THE FIRST URINE BAG WHICH CAN BE WORN LIKE A BUM BAG – WEARER COMFORT DAY AND NIGHT

Considerable improvement of mobility and thus of the patient's quality of life – discreet and simple to hide under your clothes. No change from leg to bedside bag necessary any more. Reduced risk of unintentional catheter extraction, as there is no drainage tube. Necrosis of the urethra, as caused with men by hanging urine bags, is avoided with the Belly Bag®.



Modell B 1000 CT



Modell B 1000 P

ARROW

RÜSCH

 HUDSON RCI

 Gibeck

RÜSCH CARE

MEDICAL  SERVICE

WECK

Pilling

PLEUR-EVAC

 DEKNATEL

TAUT

Teleflex serves a wide variety of customers around the world by designing, manufacturing, and distributing specialty-engineered products. Distinguished by a significant presence in healthcare, Teleflex also serves the aerospace and industrial markets with businesses that have significant market position and brand recognition.

Teleflex Medical is a global leader in disposable medical products for critical care and surgical applications and a business segment of Teleflex Incorporated (NYSE: TFX). The company also produces surgical instruments and devices, cardiac devices and other specialty products for device manufacturers. The Teleflex Medical family of brands includes the trusted names of Arrow®, Beere®, Deknatel®, Gibeck®, HudsonRCI®, KMedic®, Pilling®, Pleur-Evac®, Rüschi®, Sheridan®, SMD, SSI™, Taut™, TFX Medical OEM and Weck®.

www.ruesch.de
www.teleflex.com

Distributed by:

IRELAND TELEFLEX MEDICAL Europe, IDA Business Park, Athlone, Co Westmeath.
Phone +353 (0)9 06 46 08 00, Fax +353 (0)14 37 07 73, orders.int@teleflexmedical.com

AUSTRIA RÜSCH AUSTRIA Gesellschaft m.b.H., A-1090 Wien
Phone +43 (0)1 402 47 72, Fax +43 (0)1 402 47 72 77, rueschaustria@teleflexmedical.com

BELGIUM TELEFLEX MEDICAL BE sprl, B-1070 Brussels, Phone +32 (0)2 333 24 60
Fax +32 (0)2 332 27 40, info.be@teleflexmedical.com

CANADA RUSCH CANADA, Oakville · Ontario L6M 3E3
Phone +1 (0)905 469 8489, Fax +1 (0)905 469 8490

FRANCE TELEFLEX MEDICAL SAS, F-31460 Le Faget, Phone +33 (0)5 62 18 79 40
Fax +33 (0)5 61 83 35 84, info.teleflexmedical.france@teleflexmedical.com

GERMANY TELEFLEX MEDICAL GmbH, P.O. Box 1180, 71385 Kernen
Phone +49 (0)7151 406 0, Fax +49 (0)7151 406 150, info.de@teleflexmedical.com

INDIA TELEFLEX MEDICAL Pvt. Ltd., Chetpet, Chennai - 600 031
Phone +91 (0)44 2836 5040/-41, Fax +91 (0)44 2836 0682, mktg@teleflex-ind.com

ITALY TELEFLEX MEDICAL S. r. l., 20039 Varedo (MI), Phone +39 0362 58911
Fax +39 0362 5891 888, info.it@teleflexmedical.com

JAPAN PILLING WECK (ASIA) PTE. LTD., Minato-ku, Tokyo 105-0001
Phone +81 (0)3 5404 3415, Fax +81 (0)3 5404 3401

NETHERLANDS TELEFLEX MEDICAL B.V., P.O. Box 40, 3640 AA Mijdrecht,
Phone +31 (0)297 540 703, Fax +31 (0)297 533 726 / 540 166, info.nl@teleflexmedical.com

SINGAPORE PILLING WECK (ASIA) PTE. LTD., Lavender Place
Singapore 338750, Phone +65 (0)6294 1019, Fax +65 (0)6299 6132

SPAIN TELEFLEX MEDICAL S.A., 03530-La Nucía (Alicante)
Phone +34 (0)96 687 51 01/-02/-03, Fax +34 (0)96 587 35 69, info.es@teleflexmedical.com

SWITZERLAND RÜSCH (SCHWEIZ) AG, CH-3123 Belp
Phone +41 (0)31 819 73 90, Fax +41 (0)31 819 73 87, info@ruesch-schweiz.ch

UNITED KINGDOM TELEFLEX MEDICAL, High Wycombe, Bucks. HP12 3ST
Phone +44 (0)1494 53 27 61, Fax +44 (0)1494 52 46 50, info.uk@teleflexmedical.com

USA TELEFLEX MEDICAL, Durham, NC 27709
Phone +1 (800) 234 9325, Fax +1 (800) 932 5329, csrusch@teleflexmedical.com